

Pharmacy Plus Bluetooth Enabled Refrigerator

158L Pharmacy Plus Undercounter Fridge with Solid Door – Integrated Drawers


Product code PPSR158BT-DWP

SKU code 444411387



Key features

- ✓ Temperature range of +2°C to +8°C
- ✓ Audible and visual temperature alarms
- ✓ Solid and glass door options
- ✓ Integrated lock, supplied with 2 keys for added security

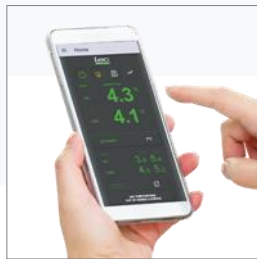

Bluetooth Connected Controller

Our controller is Bluetooth enabled and allows for simple, contactless monitoring of temperature and alarm history, via the official Lec Medical app. Battery back-up, provides continuous temperature monitoring in the event of power outage.

Lec Medical Smartphone App

Contactless local monitoring of air and load temperatures via smart phone app. Downloadable temperature report and alarm history.

Intelligent Fan Management

Smart fan optimisation aides in temperature stability following door openings, minimising fan activity when door is open to reduce the rate of temperature increase.

Dual Probe Temperature Monitoring

Our dual probe system offers greater temperature accuracy within your appliance, simultaneously measuring both air and product temperatures.

Integrated Drawers

These slide out drawers include dividers to assist with stock rotation and keep contents organised

SKU Code	444411387
Product Code	PPSR158BT-DWP
Gross Capacity (litres)	158L
Net Capacity (litres)	145L
Unit External Width (mm)	595mm
Unit External Depth (mm)	600mm
Unit External Height (mm)	850mm
Product Weight (kg)	39.86kg
Packed Width (mm)	640mm
Packed Depth (mm)	660mm
Packed Height (mm)	910mm
Packed Weight (kg)	44.91kg
Country of Origin	UK
Internal Dimensions (mm) (WxDxH)	495mm x 480mm x 720mm
Temperature Range	+2°C to +8°C
Colour	White/green
External Finish	Painted Steel
Internal Finish	HIPS*
Type of plastic	SRL800 (HIPS/ABS mix)
Door Finish	Painted Steel / HIPS*
Door Type	Solid
Door Style	Hinged, Reversible
Door Liner	Flat
Number of Door Trays	N/A
Number of Door Tray Positions	N/A
Door Open Depth (mm)	1150mm
Replacement Door Seal	Yes
Open Door Sensor	Yes
Lock	Yes, 2 Keys
Adjustable Feet (number & position)	4 All adjustable
Porthole	Yes
Porthole Size	16mm
Porthole Location	Rear Bottom (Right)
Temperature Display on Panel	Yes
Temperature Alarm	Yes
Remote Alarm Contacts	Yes
Audible Alarms	Yes
Visual Alarms	Yes
Number of Shelves/Drawers	4
Number of Shelf Positions	3 Fixed
Shelf Size (mm) (WxDxH)	N/A
Shelf Size 2 (mm) (WxDxH)	N/A
Shelf Material	N/A
Drawer Size (mm) (WxDxH)	419mm x 374mm x 60mm
Drawer Size 2 (mm) (WxDxH)	N/A
Drawer Material	ABS
Number of Baskets	1
Basket Material	Polythene Coated Wire
Basket Size (mm) (WxDxH)	485mm x 255mm x 171mm
Interior Light	Yes
Light Type	LED
Light Activation	Door Switch / App
Fan	Yes - Tangential
Spark Free	No
Defrost setting	Auto (Off Cycle)
Controller Type	Digital BT, External
Battery Back Up/Hours	Yes/48hrs
Bluetooth	Yes
Number of Probes	2
Probe Type	Air & Load (Silicone Oil)
Ambient Temperature Range	10°C to 32°C
Refrigerant Type	R600a/28g
UK 3 Pin Plug	Yes
Other Plug Types Available	No
Voltage Supply	220-240V
Frequency (Hz)	50Hz
Max Power (W)	101W
Max Current Rating (A)	0.78A
Energy Consumption (kWh/24h)	0.69kWh/24h
Energy Consumption AE (kWh)	219kWh
Noise Level (dB)	40.1dB
Insulation Thickness (mm)	45mm
Insulation Type	High Density Polyurethane Foam
Blowing Agent	Cyclopentane
Anti-Microbial Handle	Yes