

Pharmacy Plus Bluetooth Enabled Refrigerator

310L Pharmacy Plus Upright Fridge with Solid Door

Product code PPSR310BT

SKU code 444411313




Key features

- ✓ Temperature range of +2°C to +8°C
- ✓ Audible and visual temperature alarms
- ✓ Solid and glass door options
- ✓ Integrated lock, supplied with 2 keys for added security



Bluetooth Connected Controller

Our controller is Bluetooth enabled and allows for simple, contactless monitoring of temperature and alarm history, via the official Lec Medical app. Battery back-up, provides continuous temperature monitoring in the event of power outage.

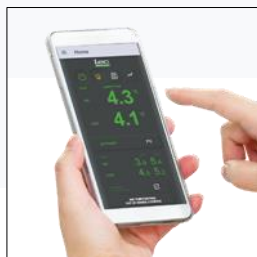

Lec Medical Smartphone App

Contactless local monitoring of air and load temperatures via smart phone app. Downloadable temperature report and alarm history.



Intelligent Fan Management

Smart fan optimisation aides in temperature stability following door openings, minimising fan activity when door is open to reduce the rate of temperature increase.



Dual Probe Temperature Monitoring

Our dual probe system offers greater temperature accuracy within your appliance, simultaneously measuring both air and product temperatures.



Flexible Internal Storage

Multiple shelf locations provides flexible storage solutions.

Product Information

Brand	Lec Medical
Category	Pharmacy
Sub-Brand	Pharmacy Plus
SKU Code	444411313
EAN Code	5052263113138
GDMA Code	PPSR310BT
Description	Free Standing Upright

Product Features

Gross Capacity (litres)	310L
Net Capacity (litres)	305L
Temperature Range	+2°C to +8°C
Open Door Sensor	Yes
Key Lock	Yes, 2 Keys
Digital Lock	No
Porthole	Yes
Porthole Size (Ø Diameter)	16mm
Porthole Location	Rear Bottom (Right)
Controller Type	Digital BT, External
Temperature Display on Panel	Yes
Temperature Alarm	Yes
Remote Alarm Contacts	Yes
Audible Alarms	Yes
Visual Alarms	Yes
Interior Light	Yes
Light Type	LED
Light Activation	Door Switch / Controller
Fan - Type	Yes - Tangential
Spark Free	No
Defrost setting	Auto (Off Cycle)
Battery Back Up/Hours	Yes/48hrs
Bluetooth	Yes
NFC	Yes
WIFI	No
GSM-4G	No
Number of Probes	2
Probe Type	Air & Load
Anti Microbial Handle	Yes

Internal Specification

Number of Shelves	5
Number of Shelf Positions	9
Shelf Size 1 (mm) (WxDxH)	496mm x 405mm x 6mm
Shelf Size 2 (mm) (WxDxH)	493mm x 265mm x 6mm
Shelf Material	White Polythene Coated Wire
Number of Drawers	n/a
Number of Drawer Positions	n/a
Drawer Size 1 (mm) (WxDxH)	n/a
Drawer Size 2 (mm) (WxDxH)	n/a
Drawer Material	n/a
Number of Baskets	1
Basket Material	White Polythene Coated Wire
Basket Size (mm) (WxDxH)	485mm x 255mm x 171mm

Door Specification

Door Finish	Painted Steel / HIPS*
Door Type	Solid
Door Style	Hinged, Reversible
Door Liner	Flat
Number of Door Trays	n/a
Number of Door Tray Positions	n/a
Door Open Depth (mm)	1150mm
Replacement Door Seal	Yes

Technical Information

Noise Level (dB)	39.9dB
Colour	White/green
External Finish	Painted Steel
Internal Finish	HIPS*
Type of plastic	SRL800 (HIPS/ABS mix)
Ambient Temperature Range	16°C to 32°C
Refrigerant Type / Weight (g)	R600a/34g
CFC Free	Yes
HFC Free	Yes
HCFC Free	Yes
Ammonia Free	Yes
Insulation Thickness (mm)	45mm
Insulation Type	High Density Polyurethane Foam
Blowing Agent	Cyclopentane
Certification	CE/UKCA
ATEX Zone Classification	n/a

Power/Installation

UK 3 Pin Plug Power Cord	Yes
EU Schuko 2 pin plug Power Cord	Yes
AU Type I Plug Power Cord	No
Voltage Supply (V)	220-240V
Frequency (Hz)	50Hz
Max Power (W)	88W
Max Current Rating (A)	0.65A
Energy Consumption (kWh/24h)	1.04kWh/24h
Energy Consumption AE (kWh)	379.6kWh
Heat Output Estimated at 22°C (kJ)	3744kJ
Adjustable Feet (Number & position)	2 adjustable (front)

Weights & Dimensions

Unit External Width (mm)	595mm
Unit External Depth (mm)	600mm
Unit External Height (mm)	1560mm
Unit External Depth Inc Handle (mm)	660mm
Internal Dimensions (mm) (WxDxH)	500mm x 480mm x 1445mm
Product Weight (kg)	51.41kg
Packed Width (mm)	640mm
Packed Depth (mm)	635mm
Packed Height (mm)	1630mm
Packed Weight (kg)	57.21kg
Country of Origin	China